


MULTIX



X5 Platform **Single Column**
110 - 150

158/180 cm

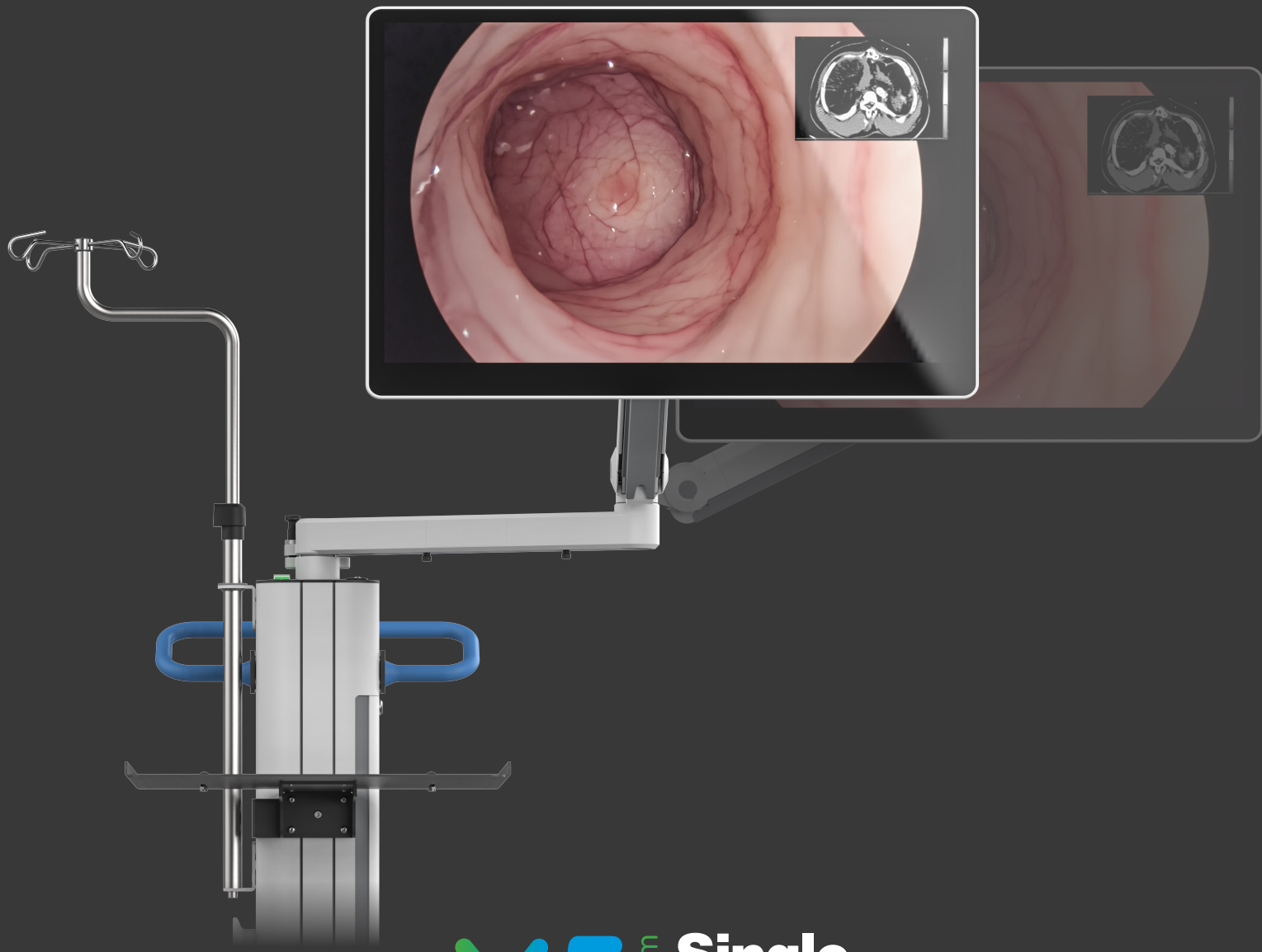

Main switch
on the top of the column

VESA Pole
for monitor up to 32"

X5 110/150

Single Column

25 kg
max load each



Tall columns for placing monitors in higher and
accommodating several devices



2/4

Height Adjustable Shelves
With an adjustable device retainer system
Inclined 2° or 8°



FOOT PRINT BASE
(W)550X(D)575 or (725 mm Long base)
TWIN WHEELS Ø 125 mm
4 with Brakes, 2 ESD


Compliant with:
MDR 2017/745 - EN 60601-1

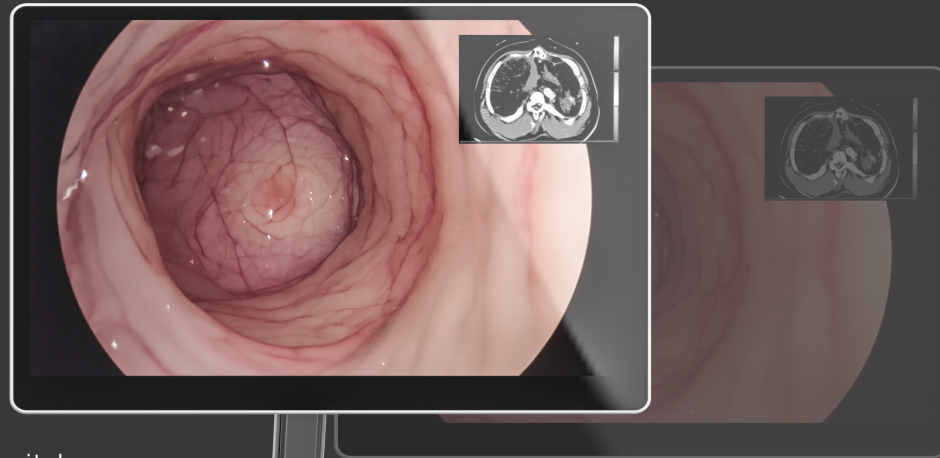
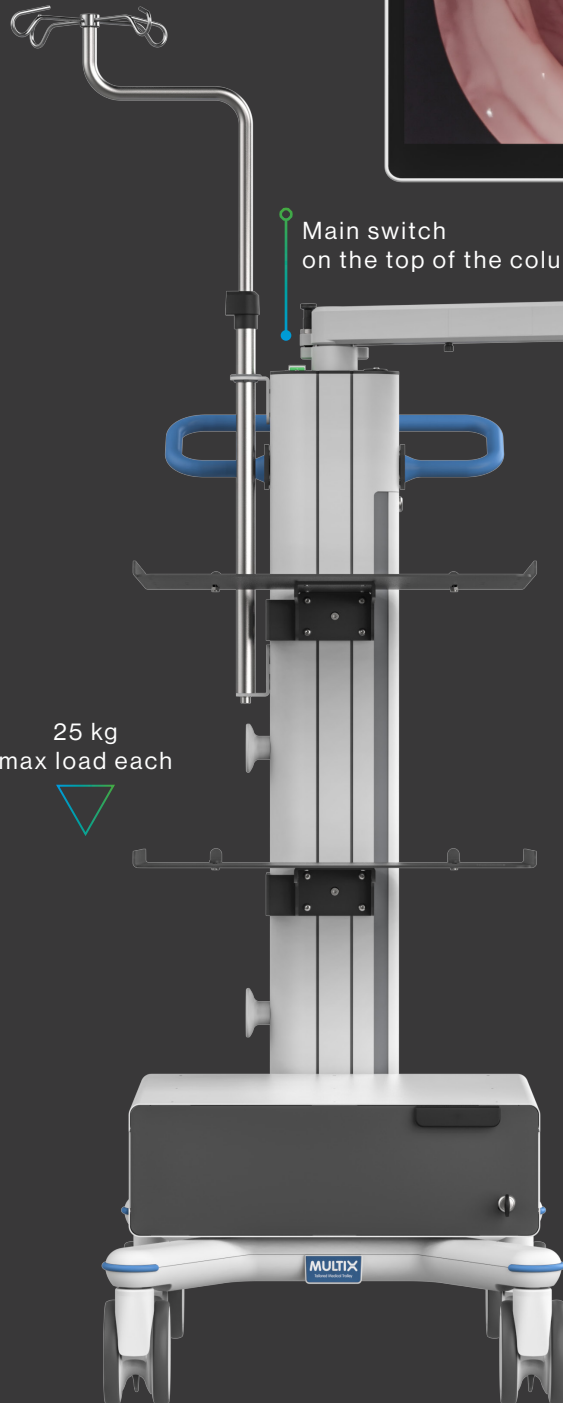



Cable management

On the X5 platform, the power strip nestles within the column, shielded by a hinged aluminium profile attached to the column itself. A flexible PVC seal facilitates cable routing from the column to the devices. Every plug has its equipotential contact and a cable clip.

It is also possible to include an isolation transformer of up to 2000W underneath the base.

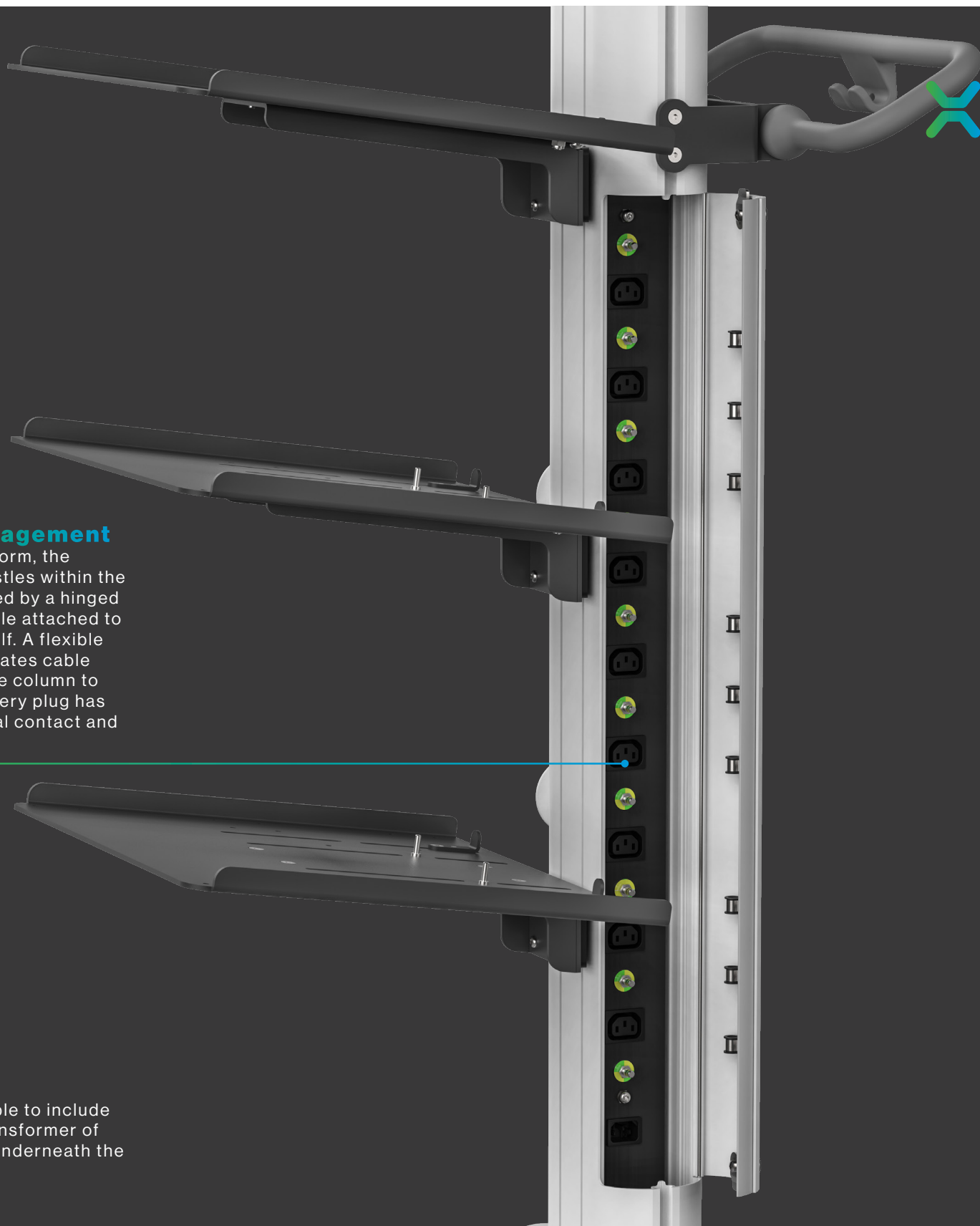

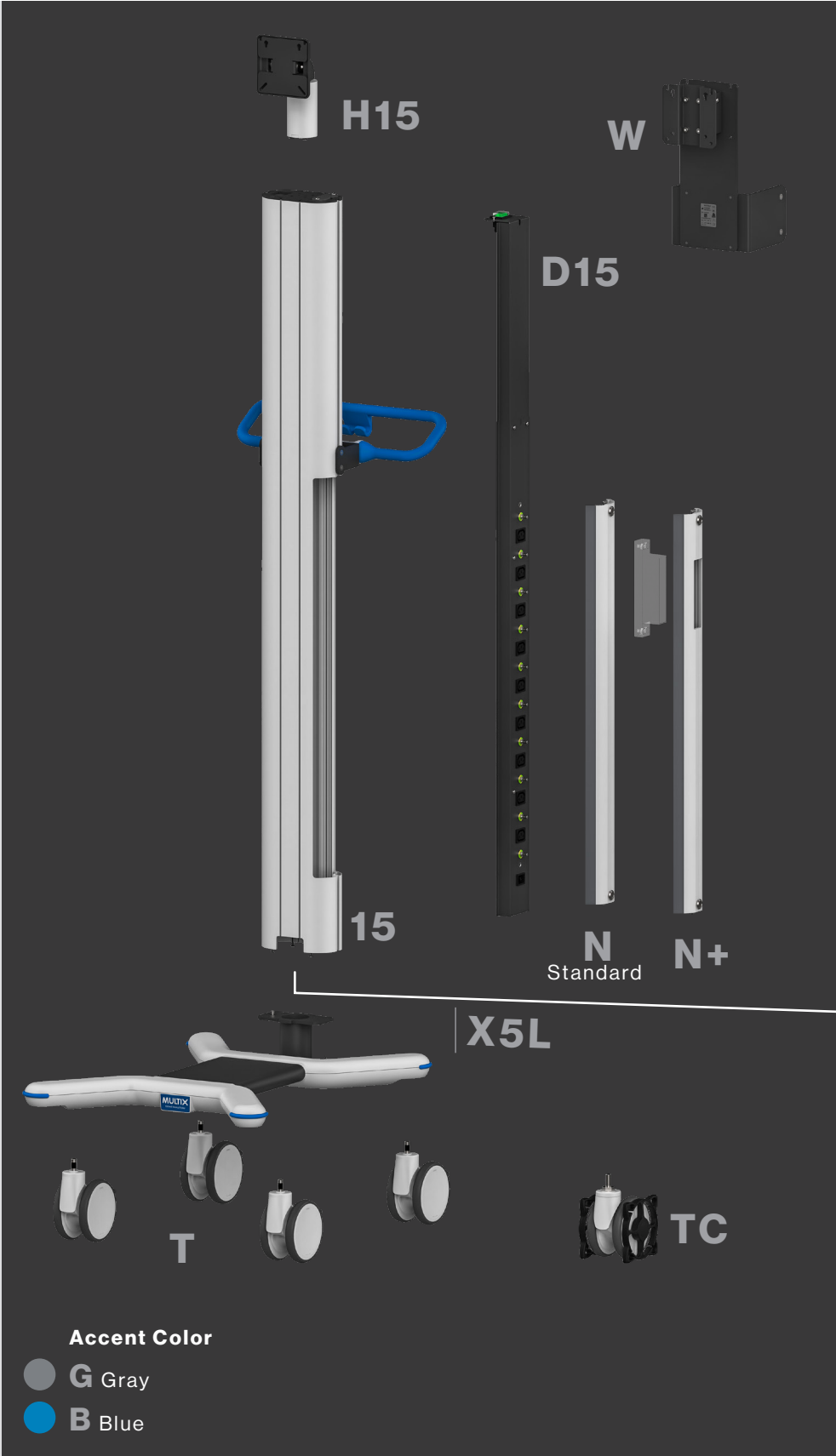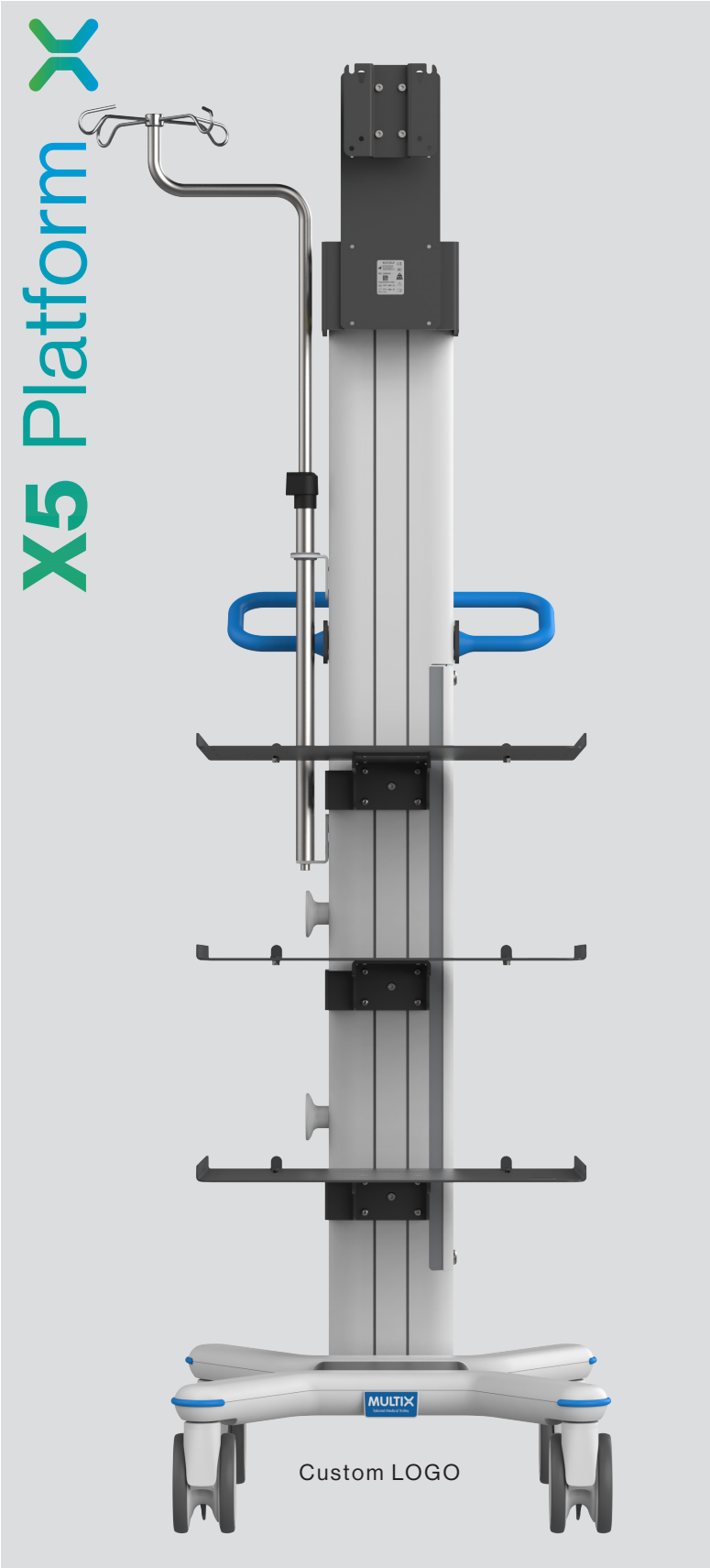

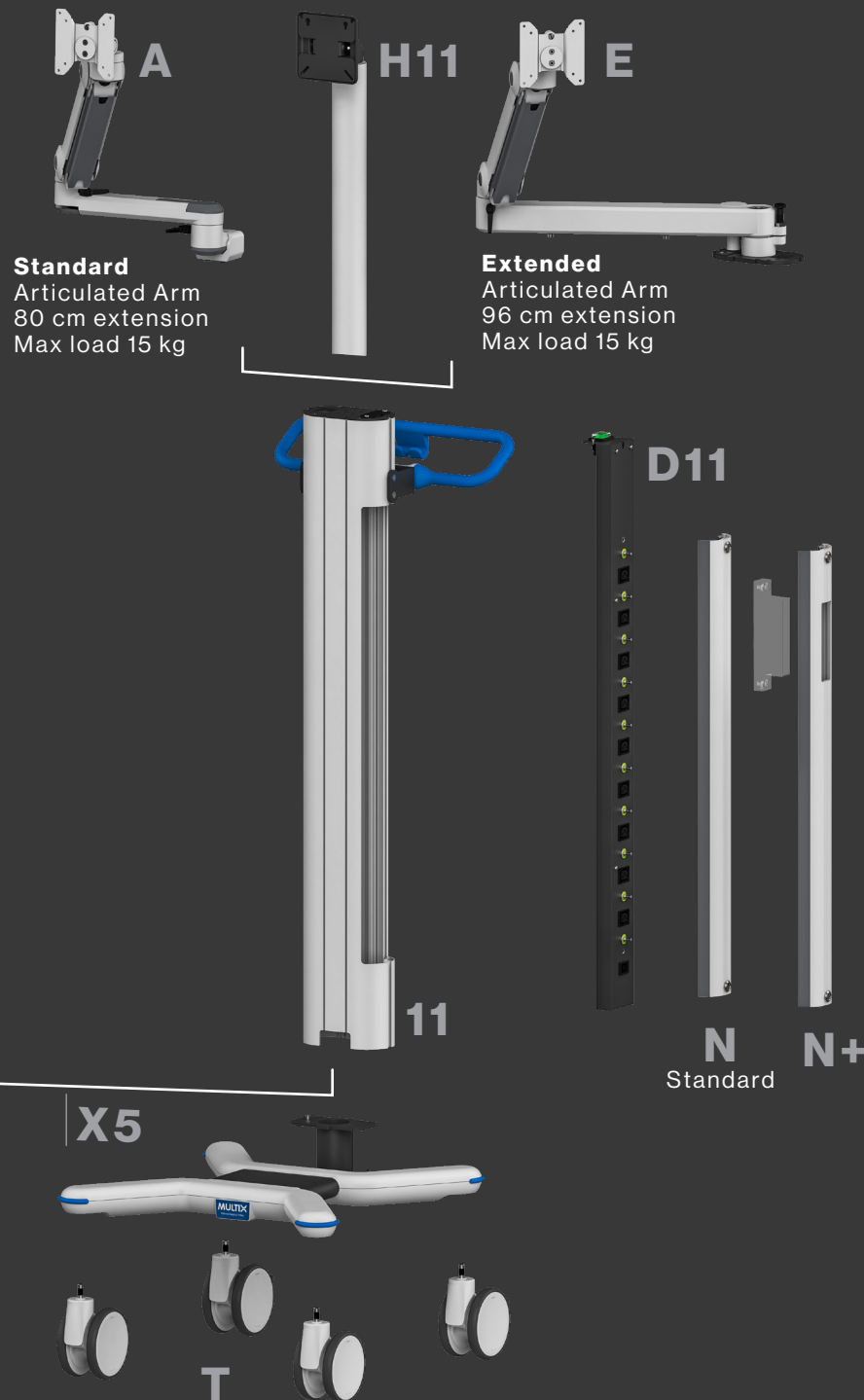


W	VESA tilting (+10- 15°) with screw system up to 55"
A	Extended articulated arm Up To 15 Kg Tilting/swivellin/Up-down
H15	VESA tilting/Swivelling Up To 16 Kg
H11	VESA tilting/Swivelling Up To 16 Kg
E	Extended articulated arm Up To 15 Kg Tilting/swivellin/Up-down
N+	Column Door With Insulator Monitor
N	Column Door Access to the plugs lockable
D 11	Multi-plug DIN
D 15	N.9 Plugs and Main Switch
11	COLUMN (H 1100mm- Total H. 1250)
15	COLUMN (H 1500mm- Total H. 1650)
X5L	FOOT PRINT (W)550x(D)725 mm
X5	FOOT PRINT (W)550x(D)575 mm
T	WHEELS Ø 125 mm 4 with Brakes, 2 ESD
TC	WHEELS Ø 125 mm 4 with Brakes, 2 ESD and Cable Guard

Platform		Color		Wheels		Column		Multiplug		Monitor Support					
X5	X5L	G	B	T	TC	11	15	D	N+	/	E	H11	H15	A	W
☐	☐	-	☐	☐	☐	-	☐	☐	☐	-	☐	☐	☐	☐	☐

X5¹⁵⁰

Example of Configuration

- > Flat shelf
- > Drawer
- > DIN rail 10×25 mm
- > IV pole for endoscopy probe support
- > Probe and camera holder
- > QWERTY keyboard and mouse support
- > Oxygen cylinder holder for 10 L

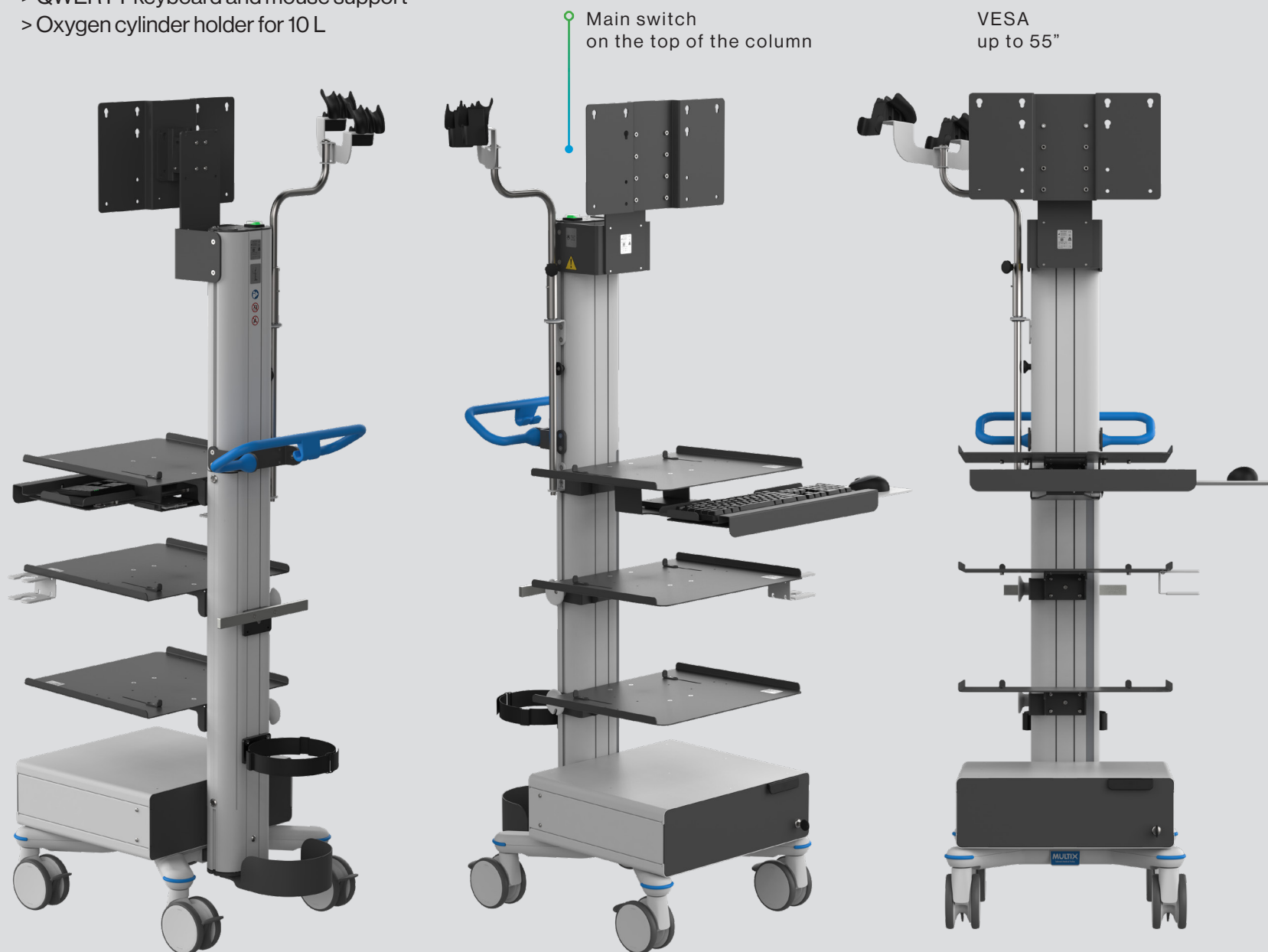



Example of Configuration


- > N. 2, 8° inclined shelf (standard or custom size)
- > Drawer
- > DIN rail 10×25 mm
- > IV pole for infusion bag
- > Oxygen cylinder holder for 5 L

Ask for all the accessories you need

- 2° inclined shelf (standard or custom size)
- 8° inclined shelf (standard or custom size)
- Drawer
- DIN rail 10×25 mm
- IV pole for infusion bag or endoscopy probe support
- Probe and camera holder
- QWERTY keyboard and mouse support
- Oxygen cylinder holder from 3 to 20 L
- Insulator Transformer up to 2000W and insulator monitor and more...



Tilting VESA
up to 32" 15kg

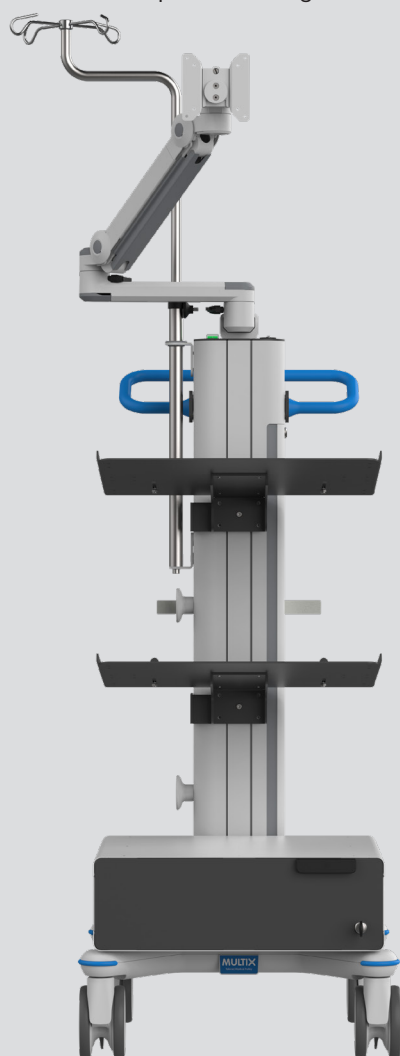



MLTCAT_X5_01(EN)



Scan the brochure




TouchPoint Group

Multix Solutions s.r.l.
phone +39 031 4100263

info@multix-trolley.com
www.multix-trolley.com

Via Padre Rovagnati, 2
22066 Mariano Comense
COMO - ITALY

Made in Italy

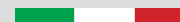

CATALOGUE **X5 TROLLEYS**