

MULTIX



X3e Platform

X3P Platform





Motorized Lifting Column

Effortless height adjustment. Our trolley features an electric lifting mechanism with a **42 cm height range**, allowing healthcare professionals to easily raise or lower the working surface to the ideal position for any user. This wide adjustment range ensures **optimal ergonomics** and helps reduce muscular strain during daily operations.

The **integrated cable management system**, fully routed inside the lifting column, keeps power and data cables protected and neatly organized. This not only improves safety and hygiene but also maintains a clean and professional appearance in clinical environments.

Thanks to its **modular configuration**, the trolley can be tailored to meet a wide variety of medical requirements through customizable compartments and accessories. This provides exceptional flexibility for securely organizing and storing essential equipment, medication, and supplies.

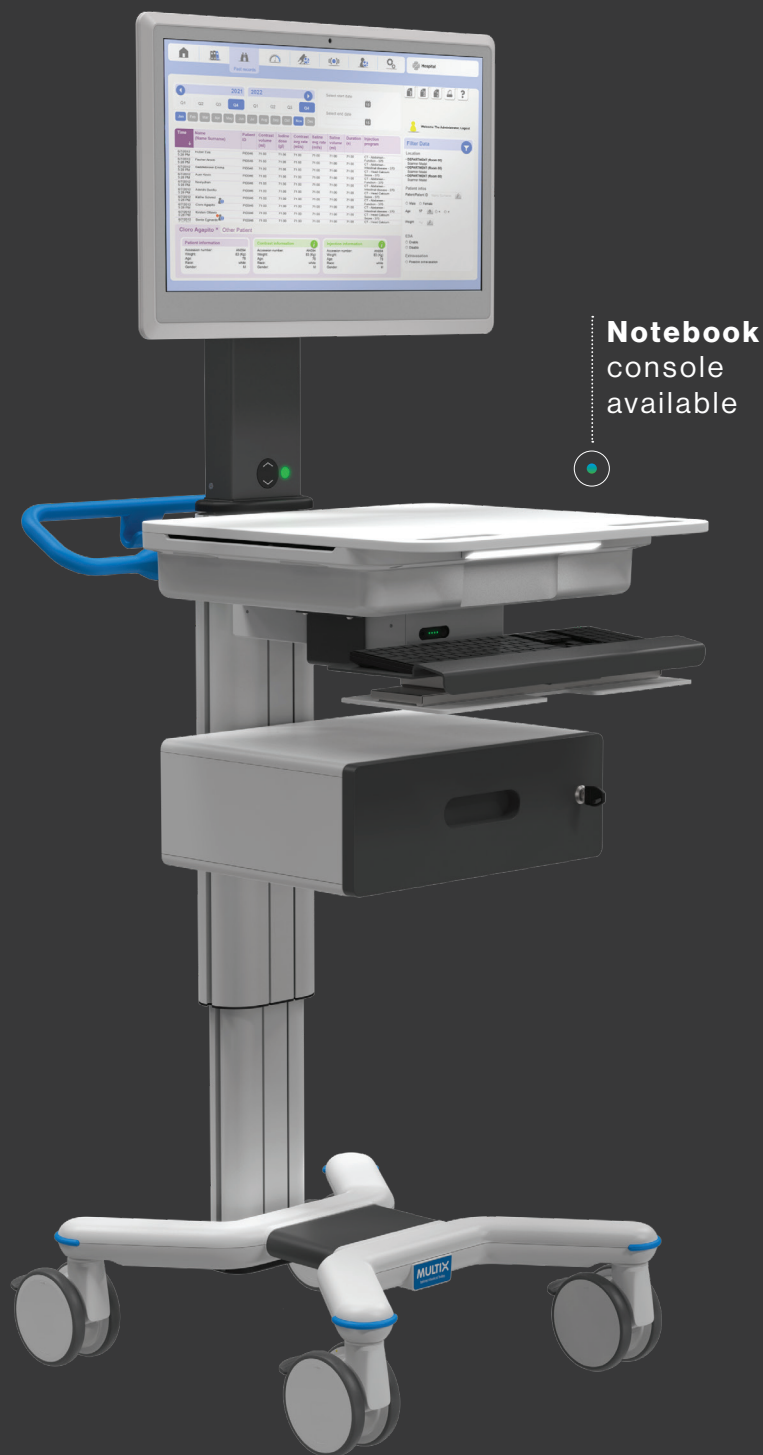


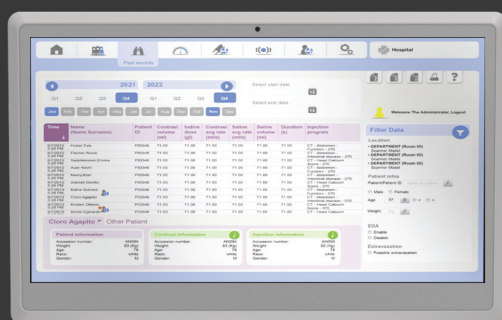


**Monitor or
Panel PC**
Up to 27"
Max Load: 12kg


**UP/Down
keyboard**

**Frontal
Handles**

UP to 42 cm

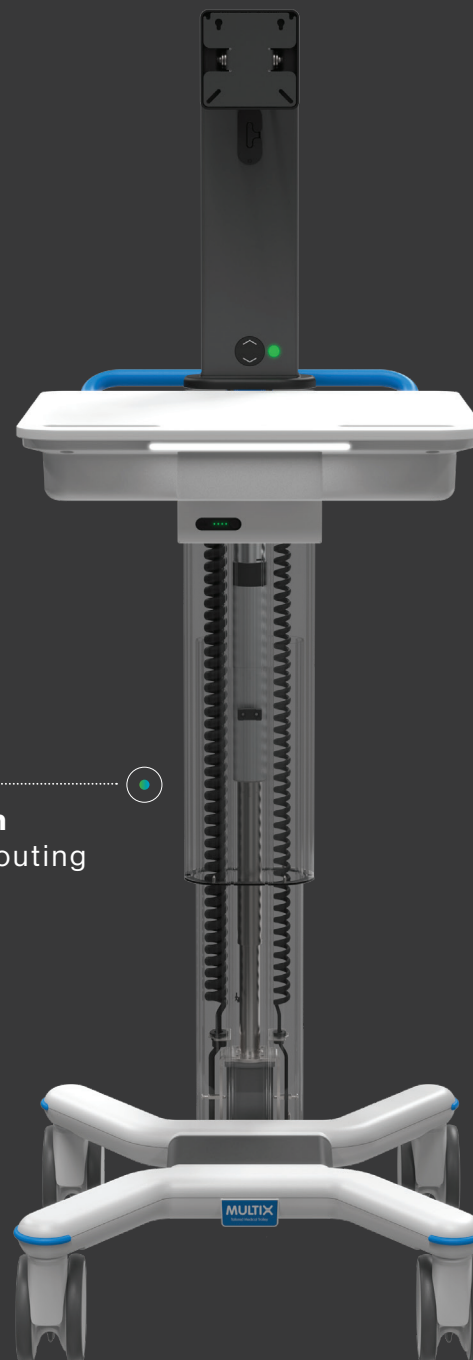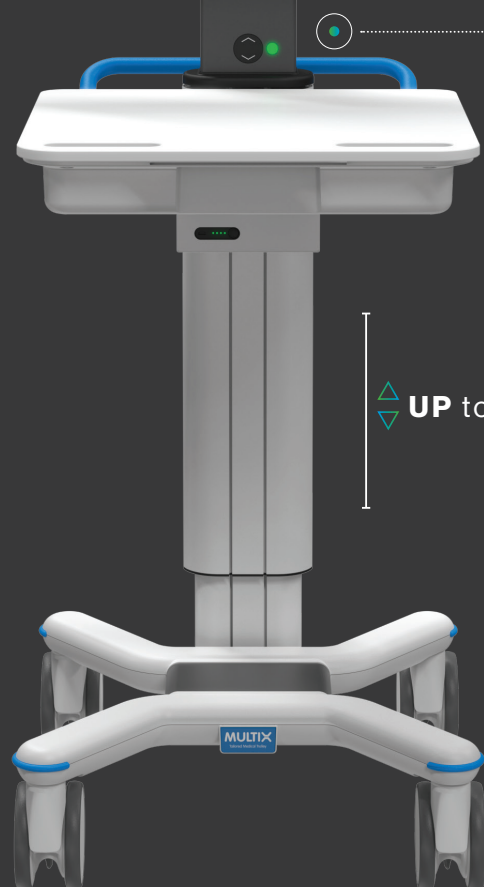

**Built-in
cable routing**

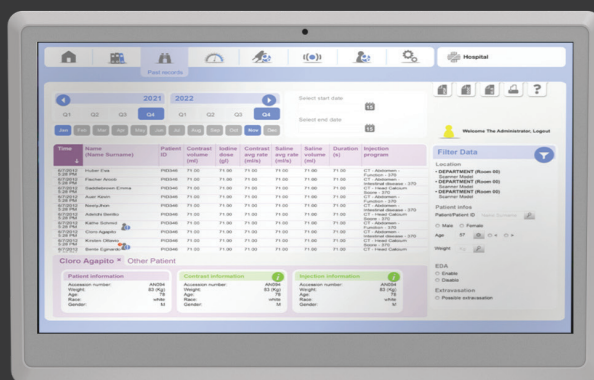

Monitor or Panel PC

Up to 27"
Max Load: 12kg

UP/Down
keyboard

Courtesy
Light

Optional
Qwerty
and MouseTray

Optional
Drawer
with standard container,
with the option to add
dividers or supplied with
ISO 40x30 trays.

Compatible with
Salto Lockers


Optional
Battery and
Display Control





Technical Box

Features five DIN connectors, a master power switch, and internal space for power supply units (PSUs) to power and manage connected equipment.

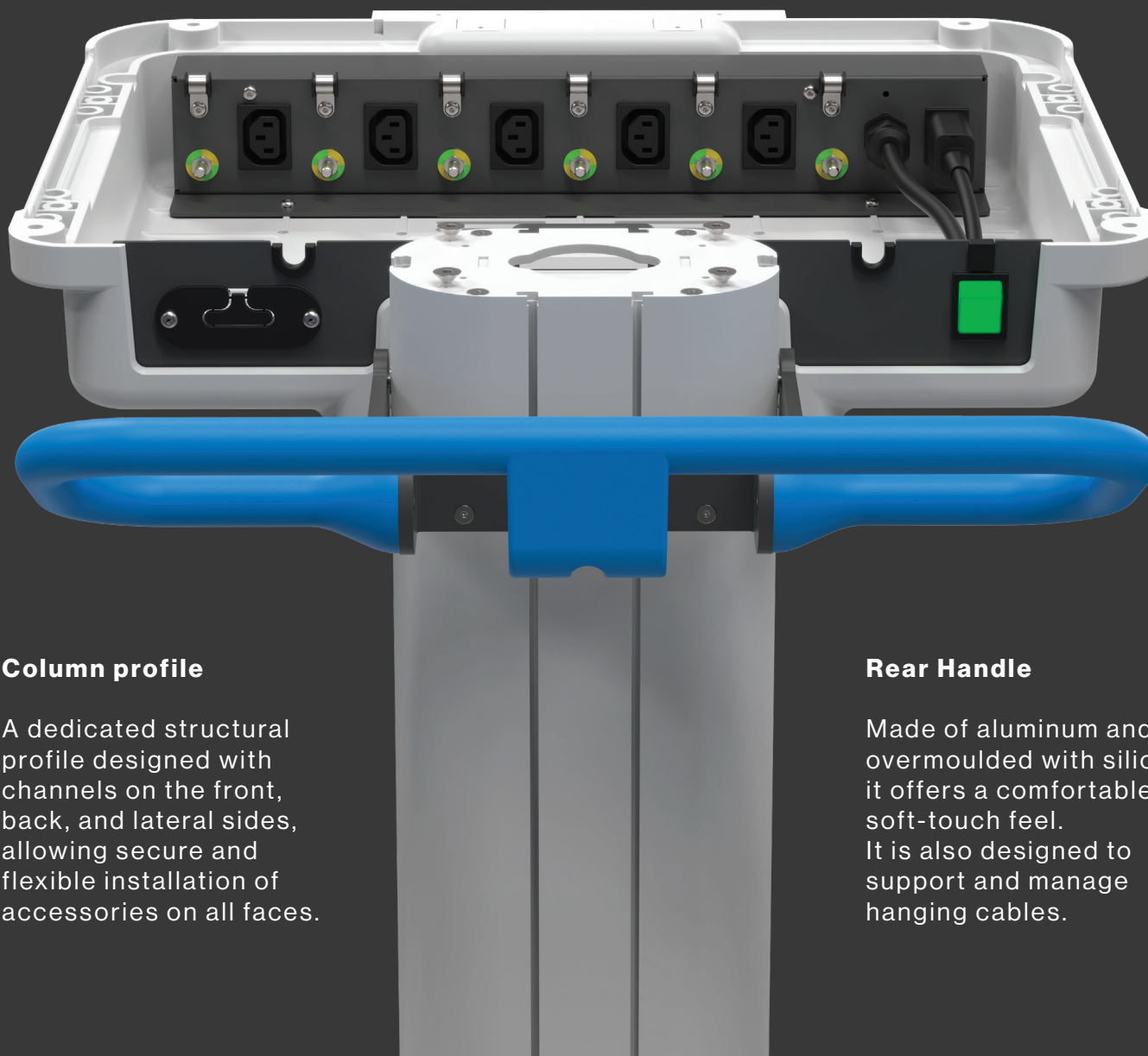

Column profile

A dedicated structural profile designed with channels on the front, back, and lateral sides, allowing secure and flexible installation of accessories on all faces.

Rear Handle

Made of aluminum and overmoulded with silicone, it offers a comfortable, soft-touch feel. It is also designed to support and manage hanging cables.




with bellow

OEM IT station

X3 Plus is equipped with **an exclusive Cable Management System**, a defining feature that sets it apart.

The bellow cable management provides an organized, clutter-free space, simplifying the integration of various healthcare devices while enhancing the overall workspace efficiency.

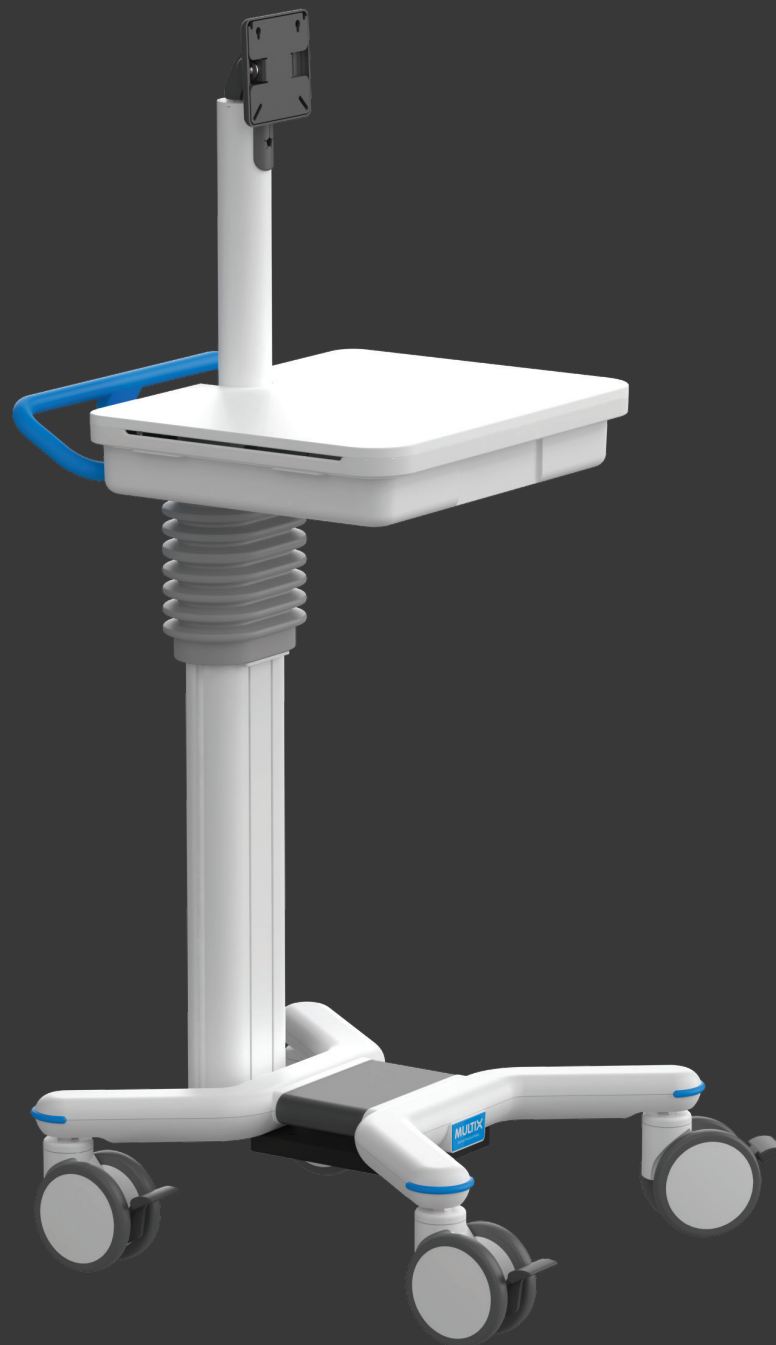

Monitor centre
from 124,5 cm to 148,5 cm




Tilting -0°+40°
Swivelling +-30°

**Monitor or
Panel PC**
Up to 27"
Max Load: 12kg

Design
The design and Ergonomy
are key to creating the
best experience for the
user.

Cleanability
Designed with smooth,
seamless surfaces,
manufactured from robust,
non-corrosive materials, and
resistant to the strongest
disinfectants.
Our trolleys meet stringent
health standards, ensuring
maximum hygiene.

Compliant with:
MDR 2017/745
EN 60601-1

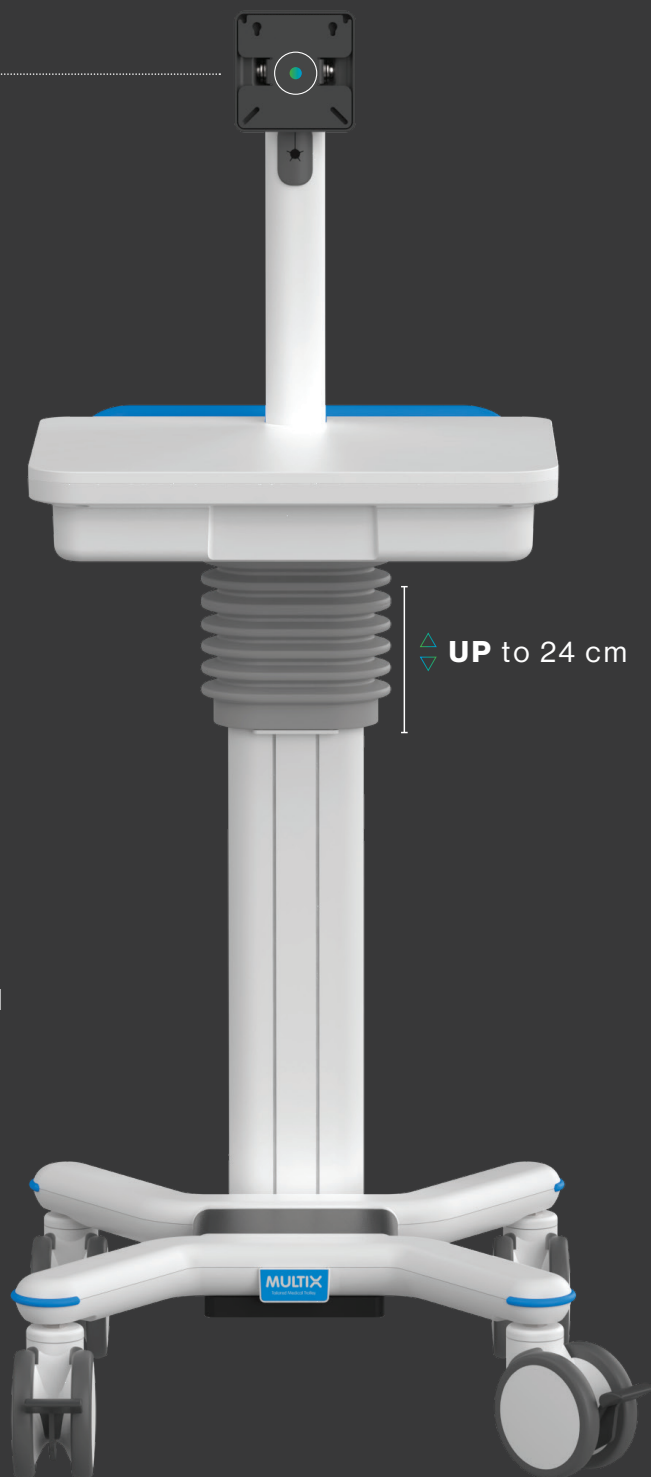

TOP PLATE
(W)513 (D)406 mm
Antimicrobial top surface on
request

INNER DIMENSIONS
(W)440 (D)280 (H) 74 mm

Work Height
from 89 cm to 114 cm
to comfortably work
standing or sitting.

**Cable
Management**
Internal cable management
systems can keep cables
contained and safe while
the trolley is moving.




FOOT PRINT BASE
(W)548 (D)578 mm
TWIN WHEELS Ø 125 mm
4 with Brakes, 2 ESD



Wide Portfolio of IT Solutions: From Notebooks to All-in-One PCS


MLTCAT X3e_X3P_01(EN)

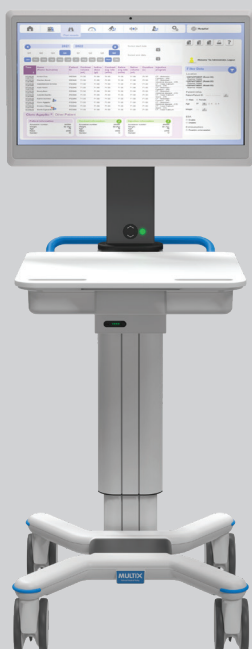


X3e
Platform




X3P
Platform

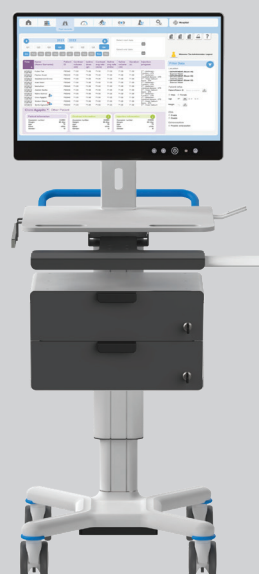

X1T
Platform

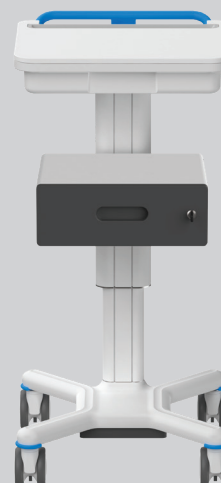

X1T
Platform

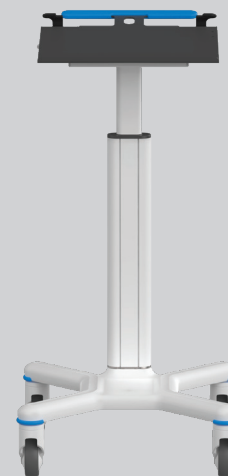

X1P
Platform



Scan the brochure



Multix Solutions s.r.l.
phone +39 031 4100263

info@multix-trolley.com
www.multix-trolley.com

Via Padre Rovagnati, 2
22066 Mariano Comense
COMO - ITALY

Made in Italy

